

## FEATURES & BENEFITS

- Single-use, disposable matrix band for use in composite restorations in dentistry.
- Can be used in all quadrants of the patient's mouth.
- No assembly or application tools required.
- Transparent band allows UV light to pass through for effective polymerisation and improves visibility of the composite
- SafeMatrix Clear Contour has the additional benefit of being pre-contoured to tooth shape.

## PRODUCT DESCRIPTION



Product name: SafeMatrix<sup>™</sup> Clear Contour

Product type: Single use, non-sterile

Material: Band: PET film, Case: ABS-Copolymer,  
Tensioning Mechanism: Nylon,  
Band Diverter Pushbutton: ABS-Copolymer,  
Adjuster Nut: ABS-Copolymer

Description: Circumferential PET film band with tightening thumb screw mechanism. A "Band Diverter" inner component with an external pushbutton helps the adjustment of the angle of the band.

Colours:

SafeMatrix Clear Contour Wide
 purple
SafeMatrix Clear Contour Narrow
 pink


Dimensions: Wide band (width): 6,0mm  
Narrow band (width): 4,5mm

Band thickness: 0,04mm

Packaging: Box of 50 pcs.

Origin: Hong Kong

For professional use only.

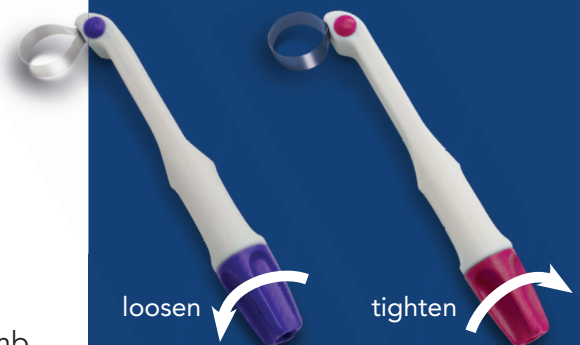

# MEDICOM® SAFEMATRIX™ CLEAR CONTOUR

Matrix Band • Medical Device, Class I

### INDUSTRIES

Dental

### CERTIFICATION & STANDARD

Meet the requirements of Regulation (EU) 2017/745 on medical devices.  
Manufacturing site ISO 13485 certified.

### STORAGE CONDITIONS

Normal conditions of conservation and storage:  
not to be exposed to moisture and sunlight.

Shelf life of the product: 5 years

### LOGISTICAL INFORMATION

#### Outer case specifications

REF	Product	Size (mm)
30065	SafeMatrix Clear Contour Narrow	280 x 243 x 300
30066	SafeMatrix Clear Contour Wide	280 x 243 x 300

#### Inner case (box) specifications

REF	Product	Size (mm)	QTY/Outer case
30065	SafeMatrix Clear Contour Narrow	130 x 115 x 135	8 boxes x 50 pcs.
30066	SafeMatrix Clear Contour Wide	130 x 115 x 135	8 boxes x 50 pcs.

