



**HEAD TO TOE  
PROTECTION**

**CATALOG**





## About us

Founded in 1988, the Medicom Group is one of the world's leading manufacturers and distributors of high-quality, single-use infection prevention and control products for the medical, dental, government and industrial markets, laboratories and research centers, animal health and health and wellness communities, and the retail trade.

## Our Vision

Medicom is a global healthcare company committed to delivering excellent products and service. We strive to partner with our customers by building strong relationships, providing optimal value, high quality products and exceptional service through continuous process improvement. We foster a work environment that brings out the best in every employee and are guided by the principles of corporate social responsibility (CSR).

## Our Mission

Medicom provides exceptional value and reliable name brand to our distributors and end users. We bring an entrepreneurial passion to everything we do to create remarkable customer experiences.

## Whom we serve

Our products are used by dentists, hygienists and other dental professionals. We offer a complete line of high-quality infection control and personal protective equipment, including a wide range of medical grade masks, gloves as well as dental specialties.

# Features & Markings



## Made in the EU

A large number of Medicom products are manufactured in the European Union in compliance with European Union hygiene and quality standards.



## Medical Device

Our products, which are classified as medical devices, meet all the requirements for being placed on the market, made available on the market and put into service in accordance with the Medical Device Regulation (MDR) or the EU Medical Device Regulation.



## Latex-free Products

This icon indicates that natural rubber latex was not used as a material in the manufacture of the product, product container, and/or packaging.

## MASKS



06

## GLOVES



14

## ESSENTIAL SUPPLIES



20

## DENTAL PRODUCTS



26

## STERILIZATION



32

# 06

## MASKS

Our clean room production 08

---

SafeMask® Medical Face Masks 09

---

Kolmi® Oxygen Respiratory Face Mask 12

---






# Our clean room production

Medicom is proud to offer a broad range of high-quality medical and FFP2 masks that provide healthcare professionals with reliable protection without compromising comfort or breathability. We manufacture our single-use masks in ISO Class 7 clean environment, which is not required by the EN 14683+AC:2019 standard, but allows us to ensure the microbiological cleanliness of our product and their ready-to-use in the operating room.

## Microbial cleanliness or “bioburden” of our masks

“Cleanliness” of the medical face mask is defined by freedom from population of viable microorganisms on a product and/or package. The EN 14683+AC:2019 standard stipulates that the medical face mask must be tested according to the EN ISO 11737-1 standard, and which provides guidance for the enumeration of the viable microorganisms on a medical device, also referred to as “bioburden”.

When tested according to EN ISO 11737-1, the bioburden of the medical mask shall be inferior to 30 CFU (Colony Forming Units) per gram of tested mask.

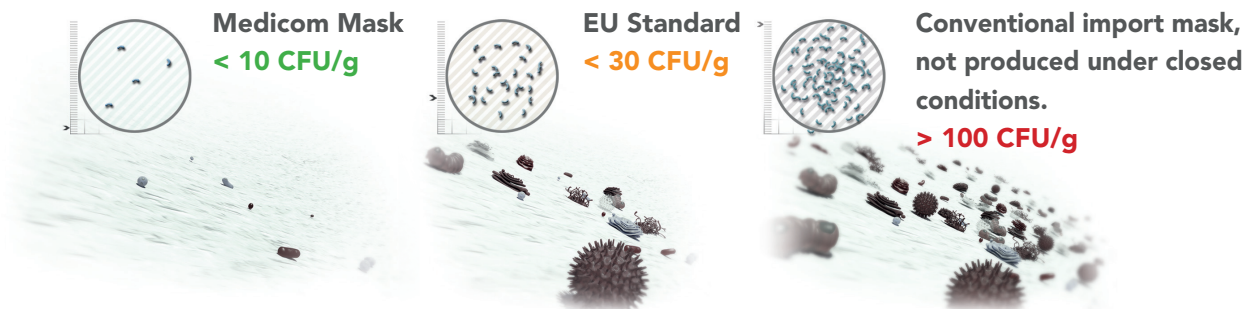

## Breathing resistance (delta P) or the differential pressure

This test determines the air permeability of the mask and its breathability. It is measured by determining the difference of pressure across the mask under specific conditions of air flow, temperature and humidity. The lower the delta P of the medical face mask, the better and easier it is for the wearer to breathe in it.

Our mask exceeds the breathability requirement of the EN 14683+AC:2019 and is one of the most breathable masks in its class, while having a high level of filtration (BFE).

# SafeMask® Medical Face Mask Typ IIR


Medicom is proud to offer a broad range of high-quality procedure masks that provide healthcare professionals with reliable protection without compromising comfort or breathability. Non-woven spunbond outer layer provides maximum protection against fluid penetration. Medicom® SafeMask® SofSkin® fog-free type IIR mask offer unsurpassed comfort and softness, making it ideal for those with sensitive skin. For those who wear eyewear, the unique “fog-free” band on the nosepiece prevents lenses and visors from fogging.

### Superior comfort

Latex-free Medicom type IIR masks offer unsurpassed comfort. SofSkin masks with a super-soft inner layer are especially recommended for users with sensitive skin.

### Splash-resistant

Medical protective masks Type IIR have excellent, liquid-resistant inner and outer nonwovens that provide a safe barrier against spray, splashes and particles.

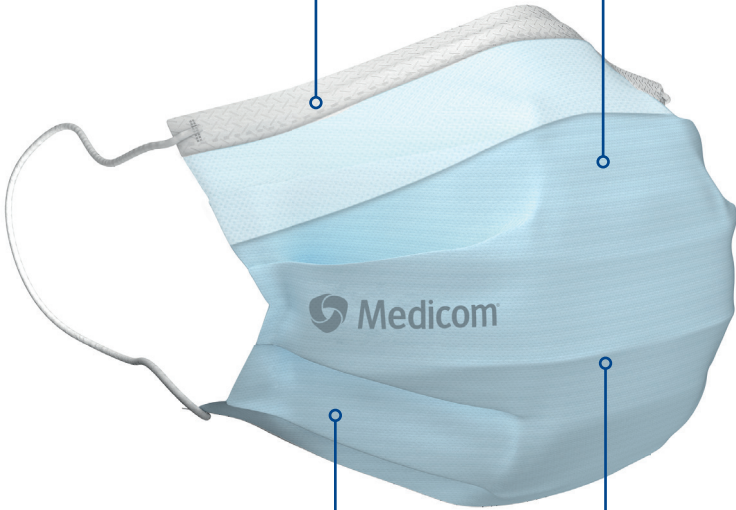

### Exceptional breathability


The breathability of the masks is unique in the market, exceeding the requirements of the standard EN 14683: 2019.

### Bacterial filtration efficiency

The high filtration efficiency (BFE) of the masks provides a safe barrier against bacterial contamination.

### The double barrier

Medical Masks Type IIR provide a physical barrier preventing other people from being exposed to bacteria and at the same time shielding the wearer from biological fluids, spray, splatter and airborne particles.

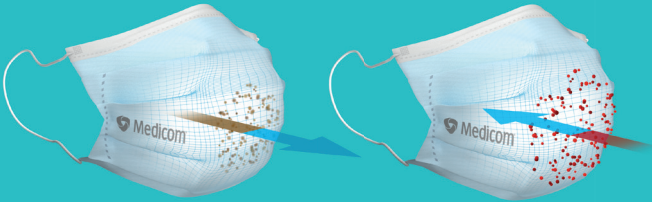







# Medical Face Masks with earloops




## SafeMask® SofSkin® fog-free with earloops

Content: 50 pcs.

	Blue	2155-FF-30
	White	2159-FF-30
	Green	2158-FF-30
	Lavender	2153-FF-30
	Pink	2156-FF-30
	Orange	2162-FF-30





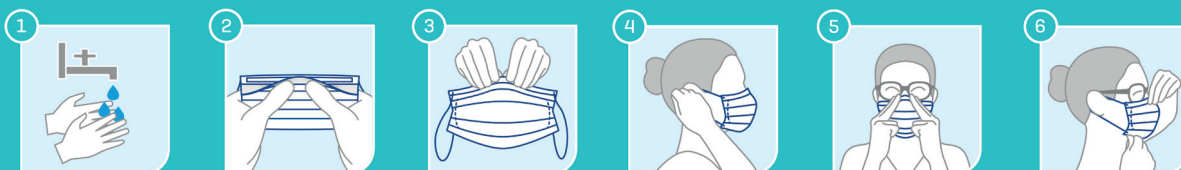



Masks are available in a "fog-free" version with a special strip on the nose bridge that prevents glasses or a visor from fogging up. To do this, lift the inner strips.



## SafeMask® Standard with earloops

Content: 50 pcs.


	Blue	2115-30
	White	2119-30
	Green	2118-30
	Lavender	2113-30
	Pink	2116-30
	Orange	2122-30



## SafeMask® SofSkin® fog-free with earloops and visor

Content: 25 pcs.

	Blue	2155-FFV-30
---	------	-------------

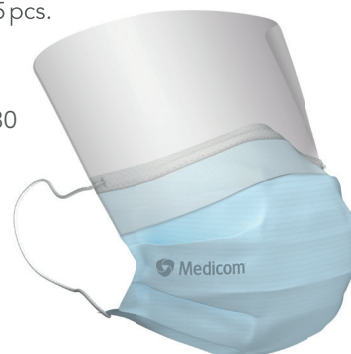

# Medical Tie-on Face Masks

In addition to a secure fit, high filtration performance, and optimal sealing, Medicom tie-on masks offer another advantage: The flexible ties allow for a customizable, pressure-free fit—ideal for extended wear. This ensures the mask remains comfortable and reliable even during demanding use.

## SafeMask® Sofskin® fog-free tie-on mask

Content: 50 pcs.


Blue 2050-BFF-30



## SafeMask® Sofskin® fog-free tie-on mask with visor

Content: 25 pcs.

Blue 2050-FFV-30

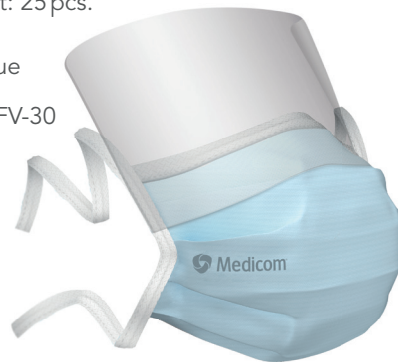

### Mask Markings



Mask with earloops



Tie-on mask



Splash resistant



SofSkin: especially for sensitive skin



with fog-free nose piece



with integrated visor

# Kolmi® Oxygen Respiratory Face Mask

High-filtration halfmask, efficiency class FFP2, type IIR. This mask effectively filters solid and liquid aerosols for reliable protection against infectious diseases (BFE 99.9%). Soft, light fleece fabric provides maximum comfort and exceptional breathability. The two ultrasonically sealed elastic head straps create a personalized fit and a tight seal around the face. Maximum 8-hour wear time.

## Individual fit

Kolmi® Oxygen FFP2 masks are available in three sizes. A special size guide helps each wearer select the correct size for optimal protection.

## Exceptionally breathable


High-quality raw materials offer exceptional breathability. Non-woven spunbond materials from leading manufacturers ensure that Kolmi® Oxygen FFP2 masks provide reliable protection and maximum comfort, even during extended wear.

## Superior comfort


Ultra-soft, latex-free inner layer makes this mask ideal for those with sensitive skin.

## High-level protection

They protect both the wearer and people in the immediate vicinity against cross-contamination. Medicom masks are also rated for efficiency and tested for overall permeability.

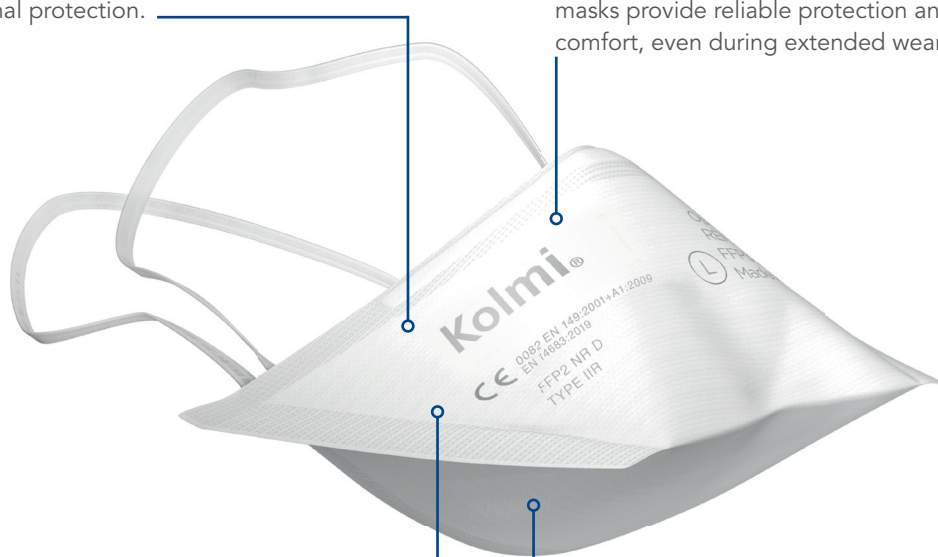

**Kolmi® Oxygen FFP2 Masks are available in three sizes to guarantee the best possible fit:**

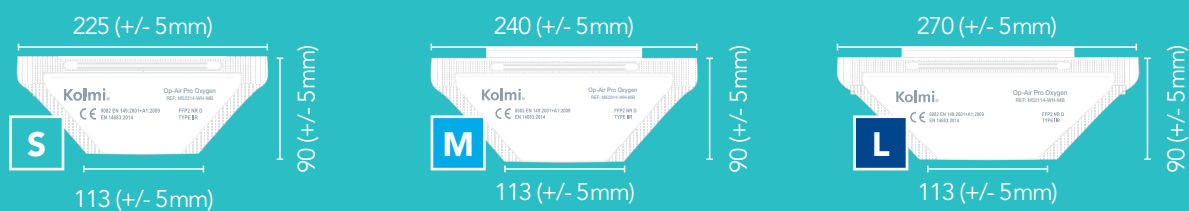



**“Well-fitting FFP2 masks reduce the risk at least into the per thousand range”**



A detailed analysis by Max Planck Institute for Dynamics and Self-Organization in Göttingen (Germany) shows that tight-fitting FFP2 masks provide 75 times better protection compared to well-fitting surgical masks and that the way a mask is worn makes a huge difference.

**S**

**Kolmi® Oxygen FFP2 Mask Small**

Content: 50 pcs.

White

**Standard**

M52214-WH-MB

**Single-packed**

M52214S-WH-MB



**M**

**Kolmi® Oxygen FFP2 Mask Medium**

Content: 50 pcs.

White

**Standard**

M52014-WH-MB

**Single-packed**

M52014S-WH-MB



**L**

**Kolmi® Oxygen FFP2 Mask Large**

Content: 50 pcs.

White

**Standard**

M52114-WH-MB

**Single-packed**

M52114S-WH-MB



# 14

## GLOVES

Latex Gloves vs. Latex-free Gloves 16

---

EU Standards 17

---

Powder free Nitrile Gloves 18

---

Powder free Latex Gloves 19

---





# Latex Gloves vs. Latex-free Gloves

Latex gloves are often preferred by dental professionals since they provide the highest level of elasticity, comfort and fit. Surgeons often prefer latex gloves because they offer greater precision with high tactile sensitivity.

Despite latex's advantages, latex allergies have become more common, so many healthcare providers are opting for latex-free glove options such as vinyl or nitrile.

Made of synthetic latex, nitrile is considered an excellent alternative for people allergic to natural latex, since these gloves best resemble the stretchy feel of latex gloves and are also highly durable and puncture resistant. Nitrile gloves provide superior protection against contamination and infections.

Another difference to consider when deciding between latex and latex-free gloves is the protection offered by each type of material. In general, latex gloves offer the best protection against bacteria and viruses, while synthetic gloves provide better chemical protection. Of the different types of latex-free gloves, nitrile gloves offer a higher degree of protection against viruses as well as protection against chemicals.

## Size table according to EN 420:

Size	Hand circumference	Hand length	Minimum glove length
6 XS	152 mm	160 mm	220 mm
7 S	178 mm	171 mm	230 mm
8 M	203 mm	182 mm	240 mm
9 L	229 mm	192 mm	250 mm
10 XL	254 mm	204 mm	260 mm
11 XXL	279 mm	215 mm	270 mm

# EU Standards

## **Medical gloves for single use – are you in safe hands?**

In Europe, a series of standards are used to secure the glove wearer and the patient from cross contamination. Below you can see a selection of the most important standards for dental.

### **EN 455-1: Medical gloves for single use**

#### **Part 1: Requirements and testing for freedom from holes**

The glove acts as a barrier. If the gloves have pinholes, this barrier and therefore the protection from cross-contamination is compromised. The test is carried out on a statistical sample. Gloves are tested during production and before shipment.

### **EN 455-2: Medical gloves for single use**

#### **Part 2: Requirements and testing for physical properties**

The standard specifies the requirements and test methods relating to the physical properties of non-reusable medical gloves, such as breaking strength, before and after accelerated aging. It is also important that the gloves are manufactured with a certain length, enabling the whole hand to be covered. The strength of the glove is important to make sure that the glove does not tear when donning the glove and using it for the intended purpose. The challenge test proves that the material will maintain its strength under certain parameters.

### **EN 455-3: Medical gloves for single use**

#### **Part 3: Requirements and testing for biological evaluation**

The glove should not cause any illness or sensitivity when wearing, and irritants should be so low as possible.

The risk of a reaction to a glove must be low. If there is a reaction, you need to be able to find out, which glove to use instead. Here the standard also helps, because it will have a list of the chemicals used in the glove.

### **EN 455-4: Medical gloves for single use**

#### **Part 4: Requirements and testing for shelf life determination**

It is important that the glove provides the expected protection. Since a glove degrades over time, it will only last for a certain period, when it is stored correctly. A glove should not be used after it is expired, as it may no longer offer the expected level protection.

### **EN ISO 374-5: 2016 Protection against micro-organisms**

The small size of a virus makes it easy to pass through the pinholes in ordinary protective gloves. Gloves that meet the EN 374-5 VIRUS Standards, however, are proven to act as an effective barrier even against small viruses. The gloves must comply with the liquidtightness test in accordance with standard EN ISO 374-2:2019.

# Powder-free Nitrile Gloves

Nitrile gloves provide a second skin feel due to their soft texture. The improved polymer technology offers a comfortable fit to provide optimal dexterity in wet and dry conditions. Textured fingertips ensure better grip control. Nitrile gloves provide high puncture resistance three times better than latex and abrasion resistance.

## SafeTouch® Advanced™ Rejuvenate

Unit weight: 4,0g (size M)  
Colour: White  
Content: 100 pcs.

6 XS 1176-A 9 L 1176-D  
7 S 1176-B 10 XL 1176-E  
8 M 1176-C



with Lanolin  
& Vitamin E

## Kolmi® Coverfeel® Touch

Unit weight: 3,5g (size M)  
Colour: Blue  
Content: 100 pcs. (XS-L)  
Content: 90 pcs. (XL)

6 XS 1140-A 9 L 1140-D  
7 S 1140-B 10 XL 1140-E  
8 M 1140-C



## SafeTouch® Advanced™ Black Lite

Unit weight: 3,5g (size M)  
Colour: Black  
Content: 100 pcs.

6 MBKN3-6XS 9 MBKN3-9L  
7 MBKN3-7S 10 MBKN3-10XL  
8 MBKN3-8M



## SafeTouch® Advanced™ Platinum

Unit weight: 3,5g (size M)  
Colour: White  
Content: 100 pcs.

6 XS 1174N-A 9 L 1174N-D  
7 S 1174N-B 10 XL 1174N-E  
8 M 1174N-C



## SafeTouch® Advanced™ Slim

Unit weight: 3,5g (size M)  
Colour: Violet blue  
Content: 100 pcs.

6 XS 1175N-A 9 L 1175N-D  
7 S 1175N-B 10 XL 1175N-E  
8 M 1175N-C



## SafeTouch® Advanced™ Vitals™




Unit weight: 3,2g (size M)  
Colour: Blue  
Content: 100 pcs.

6 XS 1193N-A 9 L 1193N-D  
7 S 1193N-B 10 XL 1193N-E  
8 M 1193N-C



# Powder-free Latex Gloves

Latex gloves are very flexible and elastic. They provide outstanding fit, comfort and dexterity and have an extreme sensitivity for good grip control. Latex offers an excellent mechanical resistance to tearing and folding and an effective barrier against biological contaminating agents.

<p><b>SafeTouch® Connect™ Rejuvenate</b></p> <p>Unit weight: 5,2g (size M) Colour: Natural (beige) Content: 100 pcs.</p> <p><b>6 XS</b> 1163-A   <b>9 L</b> 1163-D <b>7 S</b> 1163-B   <b>10 XL</b> 1163-E <b>8 M</b> 1163-C</p>  <p><b>MD</b> CLASS I</p> <p>with Lanolin &amp; Vitamin E</p>	<p><b>SafeTouch® Connect™ Soft</b></p> <p>Unit weight: 6,2g (size M) Colour: Natural (beige) Content: 100 pcs.</p> <p><b>6 XS</b> 1130-A   <b>9 L</b> 1130-D <b>7 S</b> 1130-B   <b>10 XL</b> 1130-E <b>8 M</b> 1130-C</p>  <p><b>MD</b> CLASS I</p>	<p><b>SafeTouch® Connect™ Vitals</b></p> <p>Unit weight: 5,5g (size M) Colour: Natural (beige) Content: 100 pcs.</p> <p><b>6 XS</b> 1122-A   <b>9 L</b> 1122-D <b>7 S</b> 1122-B   <b>10 XL</b> 1122-E <b>8 M</b> 1122-C</p>  <p><b>MD</b> CLASS I</p>
--	--	--

## Did you know?

Non-powdered gloves do not include cornstarch. This makes them more skin-friendly, stronger and more versatile. Powder-free gloves can be more expensive due to the process required to remove the powder coating from the glove.

However, powder-free gloves cause fewer allergic reactions because they have a lower protein and chemical content. Powder-free gloves are also cleaner to use because they do not leave residue on clothing or surfaces.



### Please note:

Latex bears a possible risk of allergies in case of prolonged and direct skin contact.

# 20

## ESSENTIAL SUPPLIES

Mouth Rinse Cups & Headrest Covers 22

---

Dental Bibs 23

---

Face Drapes 24

---














# Mouth Rinse Cups

Medicom® mouth rinse cups are made from high material density. Their 190g/m<sup>2</sup> weight ensures temperature resistance to protect hands from hot and cold. They are sturdy and stable to prevent spills. All paper cups have a 180ml capacity and fit most dispensers.

Medicom mouth rinse cups are a colourful addition to your hygiene program! Choose from eight vibrant colours and create an attractive colour scheme for your practice by coordinating with other Medicom accessories such as saliva ejectors, headrest covers and bibs.

## Mouth rinse cups

Capacity: 180ml    Height: 8,0cm  
 Thickness: 190g/m<sup>2</sup>    Content: 3.000 pcs.

 Blue	120	 Pink	125
 White	121	 Orange	126
 Green	122	 Dark blue	127
 Yellow	123	 Assorted	128
 Lavender	124		

## Printed cups

Capacity: 180ml    Height: 8,0cm  
 Thickness: 190g/m<sup>2</sup>    Content: 3.000 pcs.

 White/light beige	150
---	-----










# Headrest Covers

Hygienic headrest covers provide an additional measure of protection for the patient. The two-layer cover made of paper and foil protects against infections, prevents contamination of the headrest and provides patients with an impression of a hygienic treatment room.

## SafeBasics™ Headrest Covers

Material: Tissue + polyethylene  
 Content: 500 pieces

	25x25 cm	25x33 cm		25x25 cm	25x33 cm
 White	3014	3017	 Dark blue	3021	3025
 Blue	3015	3018	 Green	3022	3026
 Pink	3019	3023	 Orange	3028	3027
 Purple	3020	3024			















# Dental Bibs

These coloured bibs provide excellent absorption plus effective patient protection and comfort thanks to their additional paper layer and foil backing. Horizontally embossed pattern and unique water-repellent borders provide additional protection. Bibs come in a handy dispenser. Choose from a variety of bright colours to coordinate with cups and saliva ejectors.



## Dry-Back® Dental Bibs

Size: 33,0x45,5 cm  
 Material: 3-ply (2 layers of paper + 1 layer polyethylene)  
 Content: Dispenser box with 100 pcs. or 500 pcs. in carton

	100 pcs.	500 pcs.		100 pcs.	500 pcs.
 Lavender	8280-100	8280-N	 Aqua	8289-100	8289-N
 Blue	8282-100	8282-N	 Dark blue	8290-100	8290-N
 White	8283-100	8283-N	 Orange	8291-100	8291-N
 Yellow	8284-100	8284-N	 Black	8292-100	8292-N
 Green	8285-100	8285-N	 Lime	8293-100	8293-N
 Dusty rose	8286-100	8286-N	 Turquoise	8294-100	8294-N
 Peach	8288-100	8288-N			



The pack of 100 Dry-Back® dental bibs in a convenient dispenser box for easy and hygienic removal

# Face Drapes

SafeDrape™ Face drapes protect your patients during professional teeth cleaning. Whether spray mist or powder dust, the face and neck area are covered and safely protected during treatment.

The mesh fabric in the upper part allows a clear view at all times. This additional protective measure is extremely pleasant and comfortable for the patient. Female patients in particular appreciate it when the make-up is still well preserved after treatment. Available as a particularly hygienic disposable or sustainable, washable reusable product.

Today, modern dentistry means more than just professional excellence—it also stands for patient comfort and quality of service. The feeling of staying clean and protected during treatment builds trust in the practice and contributes significantly to a positive overall experience.

## SafeDrape™ Disposable Face Drapes

Size: 38,0x34,0cm  
Content: Box with 50 pieces

 Blue 8-1110B

 White 8-1110W

**MD**  
CLASS I



Available as a highly hygienic disposable option or a sustainable, washable reusable product in a convenient dispenser box.



© Sabrina Dogan, DH



The reusable SafeDrape™ face drapes are made of 100% cotton & 100% nylon and are machine washable at 40°C.








### SafeDrape™ Reusable Face Drapes

Size: 38,0x34,0cm

Material: Frottee

Content: Box with 12 pieces

	Blue	8-4020B
	White	8-4020W
	Mint green	8-4020M
	Yellow	8-4020G
	Assortment	8-4020S



# 26

## DENTAL PRODUCTS

Disposable Tips for  
the Air-Water-Syringe 28


---

Saliva Ejectors & Cotton Rolls 29

---

Disposable Matrix Bands 30


---

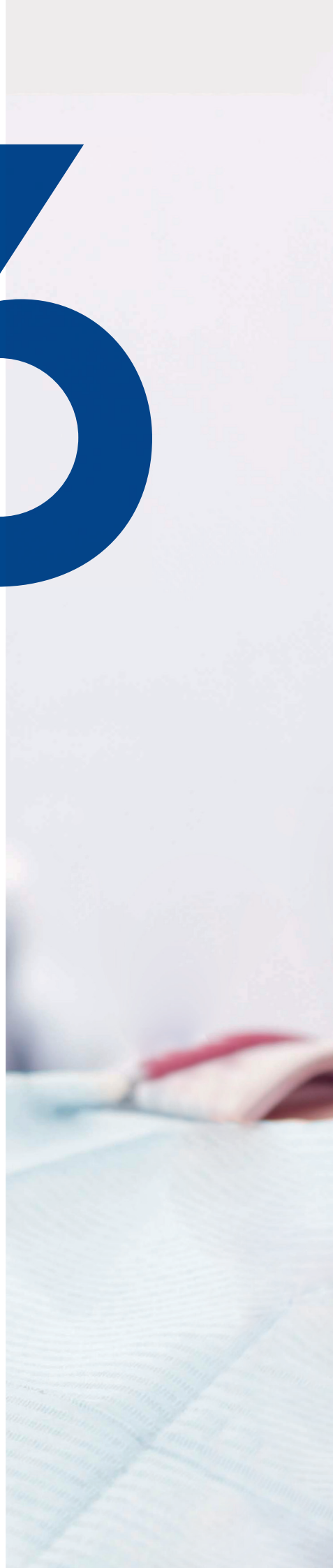



# Disposable Tips for the Air-Water-Syringe

SafeJet™ disposable tips for air and water syringes ensure infection control and prevent cross-contamination in the dental office. They provide a safer alternative to reusable syringe tips, which pose a high risk of cross-contamination. SafeJet™ ensures a reliable separation of air and water due to a unique two-channel system and produces a consistent spray. SafeJet™ fits every SafeJet™, Pro-Tip and Pro-Tip® Turbo\* converter on the market.

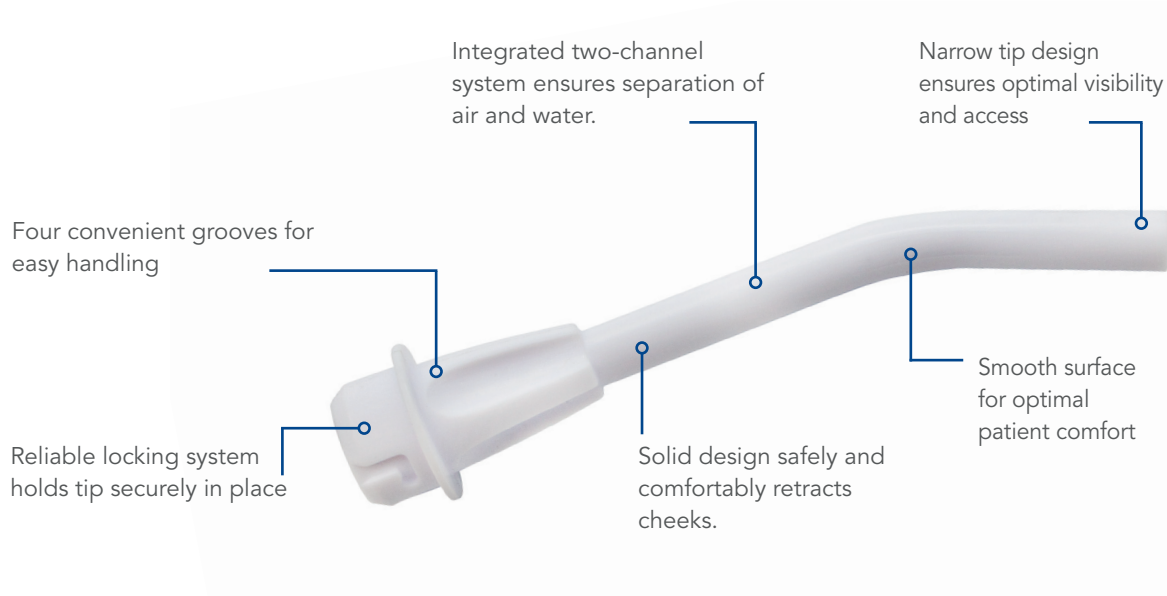

\*PRO-TIP and PRO-TIP TURBO are registered trademarks of ASTEK INNOVATIONS LIMITED. Medicom's products are not endorsed by, directly affiliated with, maintained, authorized, or sponsored by Astek.

# Saliva Ejectors








Medicom® saliva ejectors provide a continuous suction to draw saliva, blood, and debris from the mouth of the patient in order to maintain a clear operative field. The metal wire of the saliva ejector fits solidly in the tube to guarantee that the shape is perfectly maintained once bent and adapted to the patient's mouth. The rounded soft nozzle respects the sensitive gums and avoids aspiration of the flesh. The great flexibility and the absence of folds ensure an excellent aspiration.



Medicom saliva ejectors are available in lengths of 13 and 15 cm. This makes it easy to find the right saliva ejector for any situation. Choose the color that matches your practice's aesthetic!

## Saliva Ejectors – 13 cm

Content: 100 pieces

	Lavender	510-L		Yellow	508-Y
	Blue	505-B		Green	506-G
	Clear	504-C		Orange	511-O
				Pink	512-P



## Saliva Ejectors – 15 cm

Content: 100 pieces

	Lavender	710-L		Yellow	708-Y
	Blue	705-B		Green	706-G
	Clear	704-C		Orange	711-O
				Pink	712-P



# Cotton Rolls

Flexible roll maintains its shape during use. Suitable for a wide range of treatments, even on sensitive patients. Easy to remove and does not irritate the mucosa. Soft, flexible and very absorbent (absorbs 5 times its weight). Made of 100% cotton. Latex and chlorine free; non-fluorescent.

## Medicom® Cotton Rolls

Content: 300g  
Material: 100% Cotton

<b>#1</b> 0,8 x 3,8 cm	4551-300	<b>#3</b> 1,2 x 3,8 cm	4553-300
<b>#2</b> 1,0 x 3,8 cm	4552-300	<b>#4</b> 1,4 x 3,8 cm	4554-300



# Disposable Matrix Bands

SafeMatrix™ is a simple tool that makes many restorative procedures much easier for any dental team. The pre-formed matrix band ensures safe use and high patient comfort.

## Universal application

SafeMatrix is suitable for all MO/OD and MOD cavities, as well as for cusp build-ups and as a protective shield for the neighbouring tooth during preparation.

## User-friendly design

Narrow neck and rounded edges provide safe and easy access during treatment.

## Hygienic

Each matrix is individually packaged. As a single-use product, it does not require reprocessing.

## Gentle

Eliminates the risk of injuries to gingiva and cheeks often caused.

## Optimal fit


Tensioning mechanism and conical adaptation to the tooth ensure tight interproximal contact

## Compatibility

May be used with all restorative materials and with wedges and rubber dams.

## Time-saving

No need for time-consuming assembly of band and holder.

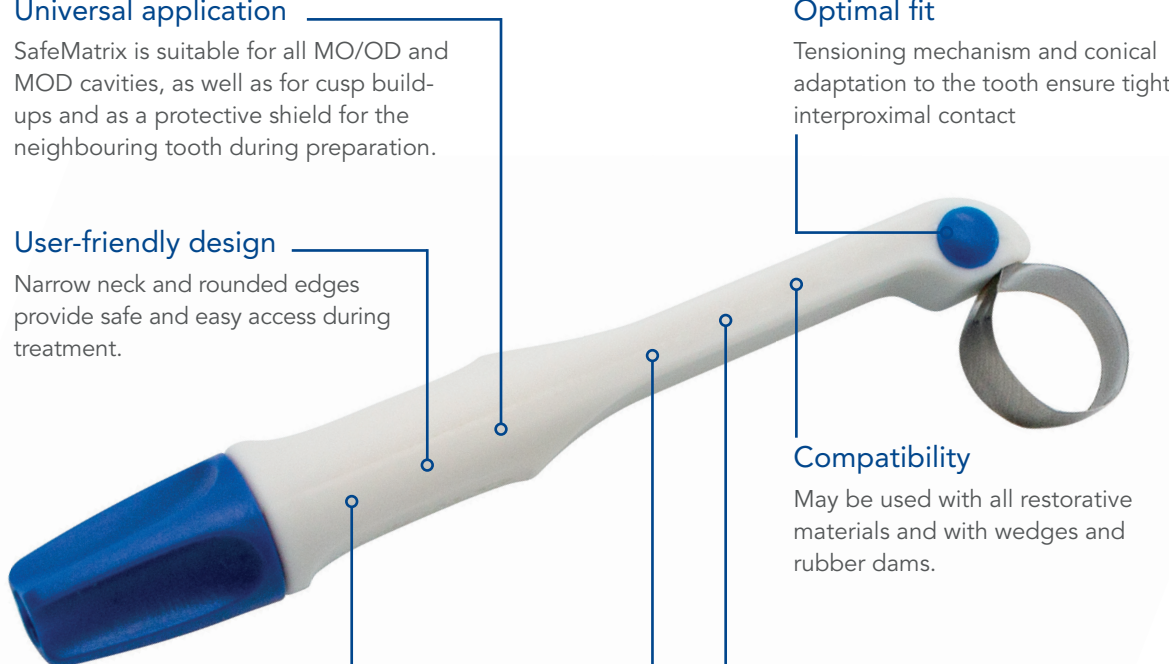

## SafeMatrix™ disposable matrix bands are available with three different types of bands

gerade

### SafeMatrix™ Standard

The popular standard version with a straight stainless steel band

konturiert


### SafeMatrix™ Contour

The anatomical curvature of the stainless steel band creates perfect contact points

konturiert

### SafeMatrix™ Clear Contour

Transparent, preformed Mylar® band allows UV light to pass through for polymerization and improves the visibility of the composite

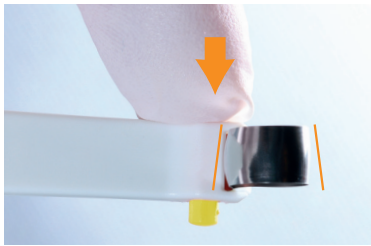

Prepare the SafeMatrix band by pushing the pin in on the side where you want to narrow the band until the pin is flush with the holder. This side must be placed towards the gingiva. The band adapts conically to the tooth anatomy.



Place the matrix band around the tooth and carefully guide it into position. Tighten the band by turning the thumbscrew clockwise. The band can be burnished if necessary. Use wedges if needed.



Create your restorative filling and cure it as usual. Remove the band by turning the thumbscrew counter-clockwise to loosen it and you're done!

© Dr. Maximilian Dobbertin

**NEW!**

### SafeMatrix™ Clear Contour

Content: 50 pcs.

**Wide** 30066

Width: 6,0mm



**Narrow** 30065

Width: 4,5mm



### SafeMatrix™ Contour

Content: 50 pcs.

**Wide** 30055

Width: 6,0mm



**Narrow** 30060

Width: 4,5mm



### SafeMatrix™ Standard

Content: 50 pcs.

**Wide** 30045

Width: 6,0mm



**Narrow** 30050

Width: 4,5mm



Extra hygienic and individually packaged in a convenient color-coded dispenser box.



# 32

## STERILIZATION

Classification 34

---

Sterilization Pouches & Reels 35

---



EO

STERILIZZAZIONE PER  
AUTOCALORE A 134°C  
PER 15 MINUTI  
NELL'ACQUA SALATA  
E STERILIZZAZIONE A  
VAPORI DI OSSIGENO  
UMIDO A 135°C PER  
15 MINUTI

STEAM

STERILIZZAZIONE PER  
AUTOCALORE A 134°C  
PER 15 MINUTI  
NELL'ACQUA SALATA  
E STERILIZZAZIONE A  
VAPORI DI OSSIGENO  
UMIDO A 135°C PER  
15 MINUTI

EO

STERILIZZAZIONE PER  
AUTOCALORE A 134°C  
PER 15 MINUTI  
NELL'ACQUA SALATA  
E STERILIZZAZIONE A  
VAPORI DI OSSIGENO  
UMIDO A 135°C PER  
15 MINUTI

Medicom

150 mm x 200 m

REF 9890

LOT

1 2 3 4 5 6 7 8 9 10 11 12

CE

MD

ISO 13485

ISO 9001

PEEL

Medicom

150 mm x 200 m

REF 9890

LOT

1 2





# Classification

Dental instruments have been classified into three categories, depending upon their risk for transmitting infection among patients and dental healthcare providers:

## Non-critical instruments

These instruments only come in contact with the intact skin. Such instruments include x-ray heads or a face-bow. These instruments have a low risk of cross-infection. Therefore, they can be cleaned and disinfected.

## Semi-critical instruments

These instruments do not penetrate oral tissues. However, they come in contact with oral mucous membrane or non-intact skin. Such instruments include dental mirrors, amalgam condensers and impression trays. It is recommended that these instruments be sterilized after each use.

## Critical instruments

These include all instruments that penetrate skin or mucous membrane, or are exposed to the bloodstream. These instruments should be sterilized after each use. Critical instruments used in dentistry include dental forceps, scalpels, bone chisels and surgical burs.

# Sterilization Pouches

Maximum efficiency! Medicom® SafeSeal® Quattro self-sealing sterilization pouches are available in a variety of sizes. Simply fill the pouch and seal it by peeling off the adhesive strip and carefully pressing down on the flap.

## SafeSeal® Quattro – Flat Sterilization Pouches

Content: 200 pcs.

<b>57 x 102 mm</b>	88015	<b>133 x 254 mm</b>	88025
<b>70 x 229 mm</b>	88000	<b>191 x 330 mm</b>	88030
<b>89 x 133 mm</b>	88005	<b>254 x 356 mm</b>	88035
<b>89 x 229 mm</b>	88010	<b>305 x 432 mm</b>	88040



Original size: 88015 (57 x 102 mm)

# Sterilization Reels

Medicom® SafeSeal® Quattro sterilization reels are available in a variety of widths that can be cut to the required length. Simply seal one end of the pouch, load the items to be sterilized and seal the other end. A sealing device is required for sealing.

## SafeSeal® Quattro – Flat Sterilization Reels

Length: 200 m

<b>50 mm</b>	9850
<b>75 mm</b>	9860
<b>100 mm</b>	9870
<b>150 mm</b>	9880
<b>200 mm</b>	9890
<b>250 mm</b>	98100
<b>300 mm</b>	98110





Medicom Europe - Dental



@medicom\_europe\_dental